YM921N

Tecnovita

by **BH**

Instrucciones de uso y garantía
Operating instructions and guarantee
Mode d'emploi et garantie
Instruções de utilização e garantia
Manuale d'uso e garanzia
Bedienungsanleitung und Garantie

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902 170 258

Thank you for purchasing our Zentromax YM921. All of our products undergo rigorous testing to ensure maximum quality. To use the unit correctly, first read the manual before using it.

INTRODUCTION


The Zentromax YM921 is based on a new technology that uses air pressure, massage rollers and vibration. It provides a soothing massage to the legs, feet and toes, reduces stress and eliminates tension in the feet. It also helps to tone and sculpt your legs. It is equipped with wheels and side handles for easy transportation.

SAFETY INSTRUCTIONS


1. Inspect the massage unit before use to ensure that it is in good condition and that the cable is not damaged. If you notice a defect or missing parts, contact the distributor of the product.
2. The unit must not be used by children or the disabled without suitable adult supervision.
3. Please place a rug under the unit to reduce the noise when the device is in operation.
4. Pregnant women should consult their doctor before using this device.
5. People suffering from a malignant tumour should consult their doctor before using this device.
6. People suffering from headaches should consult their doctor before using this device.
7. People who have recently undergone medical surgery should consult their doctor before using this device.
8. Do not use it directly on varicose veins or open wounds.
9. Do not use this device outdoors.
10. Do not handle the plug with wet hands.
11. Do not place the device in places that are too close to heat sources.
12. Do not use the device in dusty, oily or smoky environments.
13. Do not use the device for more than 30 minutes.
14. Do not use this device within 20 minutes after meals.
15. Do not use the massager in damp/humid surroundings and prevent it from getting wet. Do not touch the appliance with wet hands while it is operating. Under no circumstances can this product be used in the bath or shower.
16. Do not grip the appliance by the cable.
17. Do not expose the massager to high temperatures. Store it in a dry place.
18. Do not unplug the appliance without switching it off first.
19. To prevent the risk of fire or electrical damage, use the correct voltage.
20. In the case of a power outage unplug the device from the mains supply.

21. If you notice excessive heat or noise turn the device off and wait at least 1 hour before using it again.
22. Unplug the device from the mains if it is not going to be used.
23. Avoid knocking the device and treat it as an electrical product
24. Make sure the On/Off switch is in the Off position before plugging it into the mains supply.
25. Unplug the device before cleaning it, use a soft dry cloth to clean it. Do not use abrasive cleaning products.
26. If you notice that the massager is not operating correctly, do not attempt to repair, dismantle or modify it. Unauthorised repairs may cause injury, in addition to breaking the terms of the warranty. Contact your technical assistance service.

DESCRIPTION

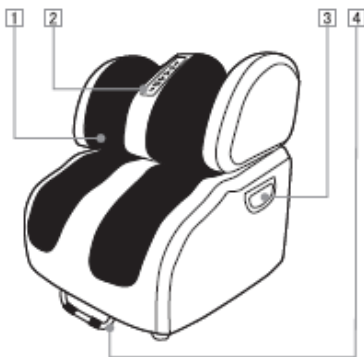
1. **Cloth cover**
2. Control panel
3. Handle
4. Front support

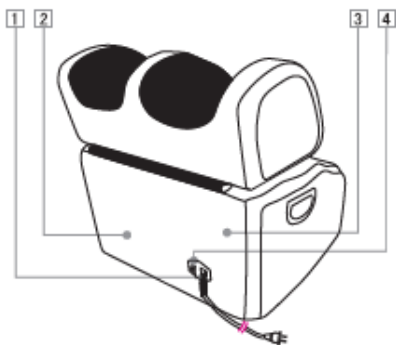1. Input
2. Rollers
3. Cover
4. On/Off switch

CONTROL PANEL

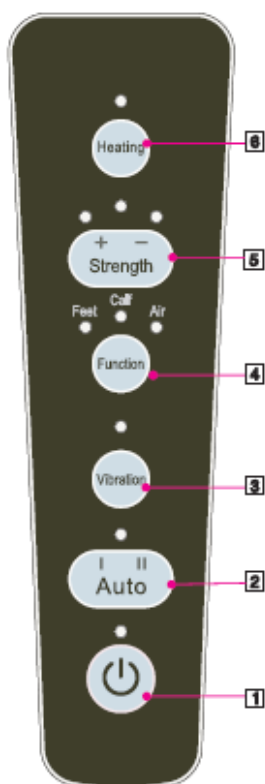6 – HEATING

To select low temperature, press the button once, the indicator light will stay on. To select high temperature, press the button again, the indicator light will stay on flashing. To switch the heat function off press the button again.

5 - STRENGTH

Press + to increase the intensity of the massage Press - to decrease the intensity of the massage

4 – FUNCTION

Press the button once to select foot massage, the FEET light will come on. Press the button again to select leg massage, the CALF light will come on. Press the button again to select air massage, the AIR light will come on.

3 – VIBRATION

Press the button to turn on the vibrating massage. This function can be used in automatic mode, during FEET or CALF mode.


2 – AUTO

Automatic modes. Press I to choose mode I, focusing on the foot area. The indicator light will remain lit. Press II to choose mode II, focusing on the leg area. The indicator light will flash.

1. – ON/OFF BUTTON

Press the button to switch the unit on. The indicator light will light up. Press the button again to switch the device off. The indicator light will switch off.


INTENSITY LEVEL DIAGRAM

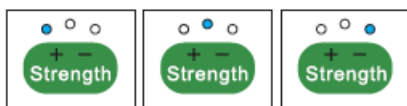1

2

3

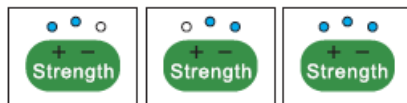4

5

6

OPERATION

1. Lift the front support before using the device.

2. Plug the cable into the mains supply

3. Press the On/Off switch at the back of the device

4. Press the On/Off button on the control panel.

5. Select the desired massage settings by using the buttons on the control panel.

For best results, 15-30 minutes continuous use is recommended.

The appliance will switch off automatically after 15 minutes of use.

NOTE: Before plugging the cable into the mains supply, make sure the on/off switch at the rear of the unit is in the 0 position.


If you do not notice any heat when the heat function is on, make sure the cable on the heat pad is not loose.

MAINTENANCE

1. Make sure the device is not connected to the mains supply before attempting to clean it.
2. Clean the device with a damp cloth (chemicals such as detergents, alcohol or solvents should not be used to clean the device).
3. Use a dry cloth to dry the device.
4. Allow the device to dry naturally, do not use other methods to dry the device quickly.
5. Do not store the device in high temperatures, do not expose the device to direct sunlight or where there are corrosive gases.
6. Keep the device away from moisture, dust or smoke.



The cloth covers can be removed for easier cleaning.

TRANSPORTATION

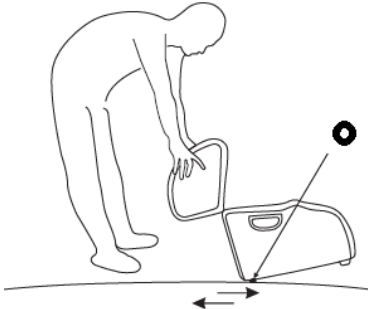TECHNICAL SPECIFICATIONS

MODEL	YM921 ZENTROMAX
POWER	120W
DIMENSIONS	540x530x530(mm)
VOLTAGE	220V~50Hz
WEIGHT	21.3kg

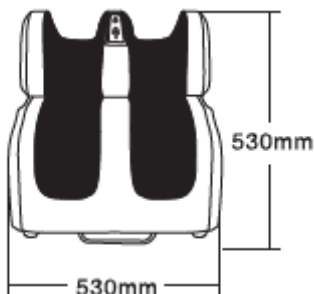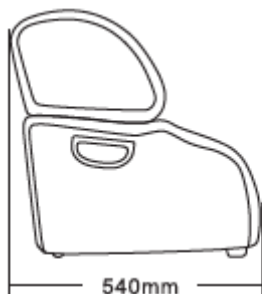Warranty

This product is guaranteed for 2 years in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25th May, 1999, regarding certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10th July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT n^o B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product.

Date of purchase.-	Purchased in.-
Purchaser.-	Dealer stamp.-

Please keep the receipt as proof of purchase.

Should you have any queries please do not hesitate to phone our customer service line:

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BH China Co., Ltd.