Tecnovita
by **BH**

Instrucciones de uso y garantía
Operating instructions and guarantee
Mode d'emploi et garantie
Instruções de utilização e garantia
Manuale d'uso e garanzia
Bedienungsanleitung und Garantie

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902 170 258

Thank you for purchasing our YR30 bt belt electrostimulator.
All of our products undergo rigorous testing to ensure maximum quality.
To use the unit correctly, first read through the whole manual before using it.

INTRODUCTION


EMS, electric muscular stimulation is a widely recognized method that has been used on the sports field for a long time. The bt belt electrostimulator is both easy and comfortable to use. This product is specially designed to help tone the abdominal muscles. It allows you to enjoy the exercise and increase muscle capacity by developing strength and improving tone effortlessly. Thanks to the design of the bt belt electrostimulator it can tone the abdominal muscles anytime, anywhere.

GENERAL SAFETY INSTRUCTIONS


1. Before using the appliance make sure it is in good working order. If you detect any faults get in touch with your local dealer.
2. The unit must not be used by children or the disabled without suitable adult supervision.
3. Pregnant women must not use it.
4. Do not use if you suffer from heart disease or other type of disease.
5. Do not use the appliance if you have cancer.
6. Do not use if you are prone to attacks, for example, epilepsy.
7. Do not use the appliance after surgical operations in which increased muscle tension may affect the recovery process
8. Do not use directly on open wounds or varicose veins.
9. Do not use if you have any implanted electrical device, such as a pacemaker.
10. Do not use if you have any metal implants.
11. Do not use near the genitals.
12. Do not use on broken or inflamed skin.
13. Do not use the appliance near your heart.
14. If you suffer from any serious illness, especially if you suspect that you may have high blood pressure or problems with blood clotting, consult your doctor before using the appliance.
15. Consult your doctor before using this product if you have any type of skin disease.
16. If you suffer from chronic pain of unknown origin, regardless of where on your body, consult your doctor before using the appliance.

17. People undergoing medical treatment should consult a doctor before using this appliance.
18. Diabetics should consult their doctor before using this appliance.
19. People with constant skin irritation resulting from prolonged stimulation by the same electrode, should consult a doctor.
20. Consult your doctor if you have any problems resulting from the stimulation treatment.
21. Stop using the appliance if you experience any discomfort during use.
22. Do not use after consuming alcohol.
23. If the appliance malfunctions or if you feel discomfort or pain, discontinue use immediately.
24. Do not use while sleeping, driving or while using some other device.
25. Keep the appliance away from heat sources and do not use close (<1 m) to medium or microwave wave devices (e.g. mobile phones) because they can cause power surges.
26. Do not use when there are high frequency surgical devices connected simultaneously.
27. Do not allow the electrodes to come into contact with metal objects, e.g. belts or collars, during operation.
28. Do not use near flammable, gaseous or explosive substances.
29. Switch the appliance off before removing it to prevent unnecessary irritation.
30. Never use the appliance without first applying plenty of electro-stimulation gel. If you do, then the stimulating effect will not be transmitted correctly and it can also damage the appliance.
31. Do not apply enriching creams beforehand, as these cause the electrodes to wear out faster or may cause unwanted current surges.
32. Do not use the unit in damp/humid surroundings and prevent it from getting wet. Do not touch the device with wet hands while it is operating. Under no circumstances must this product be used in the bath or shower.
33. Do not expose the appliance to direct sunlight or high temperatures.
34. Do not allow the appliance to come into contact with water or any other liquid. Do not immerse the appliance in liquid to clean it.
35. If you notice that the appliance is not operating correctly, do not attempt to repair, dismantle or modify it. Unauthorised repairs may cause injury, in addition to breaking the terms of the warranty. Contact your technical assistance service.

DESCRIPTION

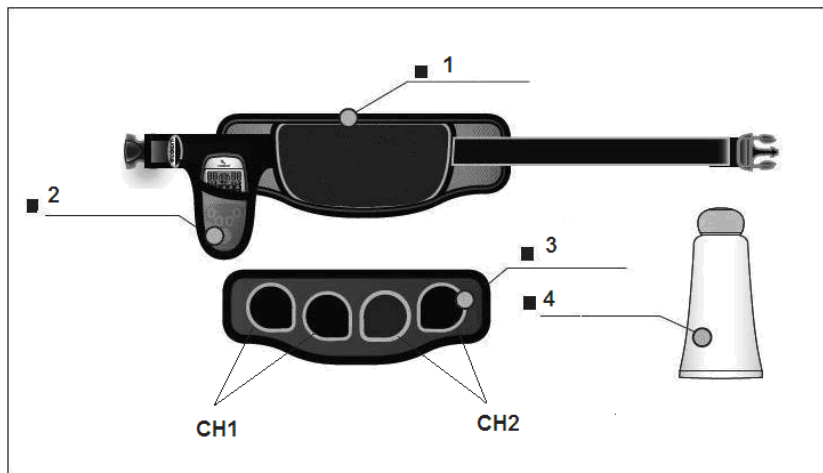
- 1 - BELT
- 2 - CONTROL UNIT
- 3 - 4 ELECTRODES
- 4 - ELECTRO-STIMULATION GEL
- CH1 - CHANNEL 1
- CH2 - CHANNEL 2

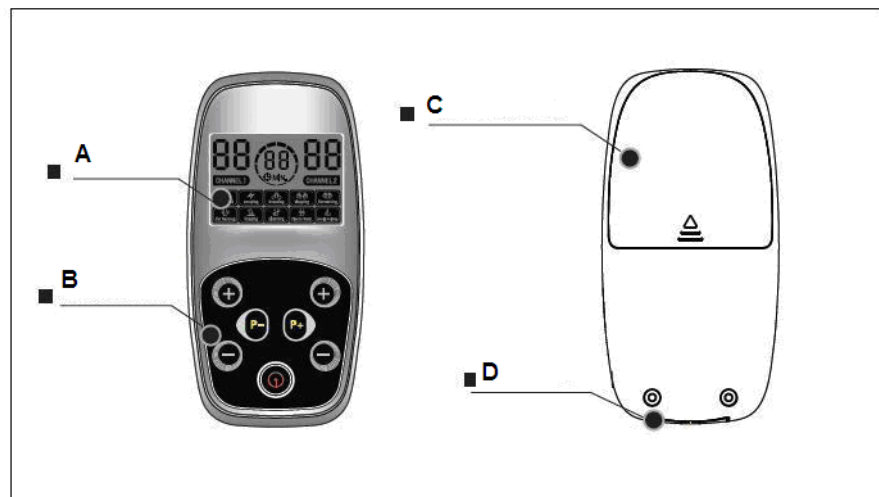- A - Digital screen
- B - Operating buttons
- C - Battery compartment cover
- D - EMS output

OPERATING INSTRUCTIONS.-

CONTROL UNIT

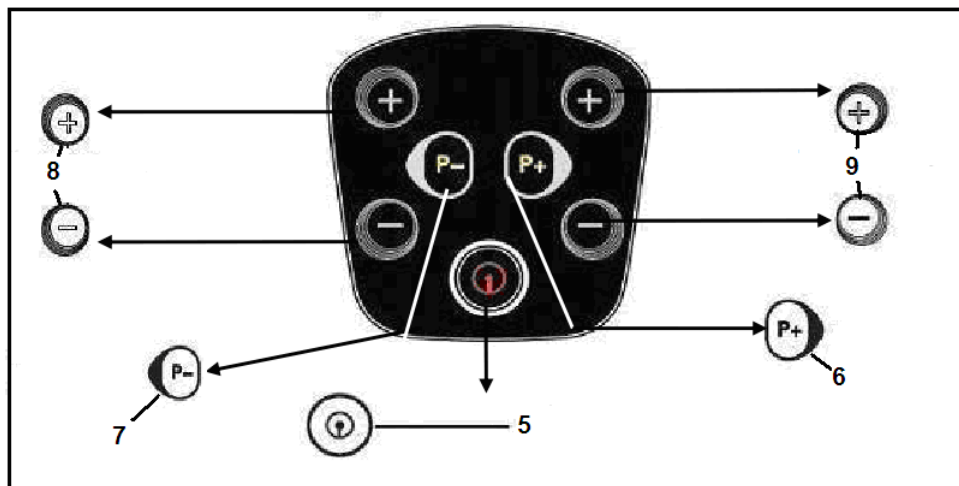5 - On/Off button


6, 7 – Mode selection (10 modes)

8 – Channel 1 intensity setting

9 – Channel 2 intensity setting

	Press to switch appliance on Press and hold to switch appliance off
	(Left side) Channel 1 setting. Press +/- to increase or lower intensity.
	(Right side) Channel 2 setting. Press +/- to increase or lower intensity.
	Select backward mode.
	Select forward mode.

DIGITAL SCREEN

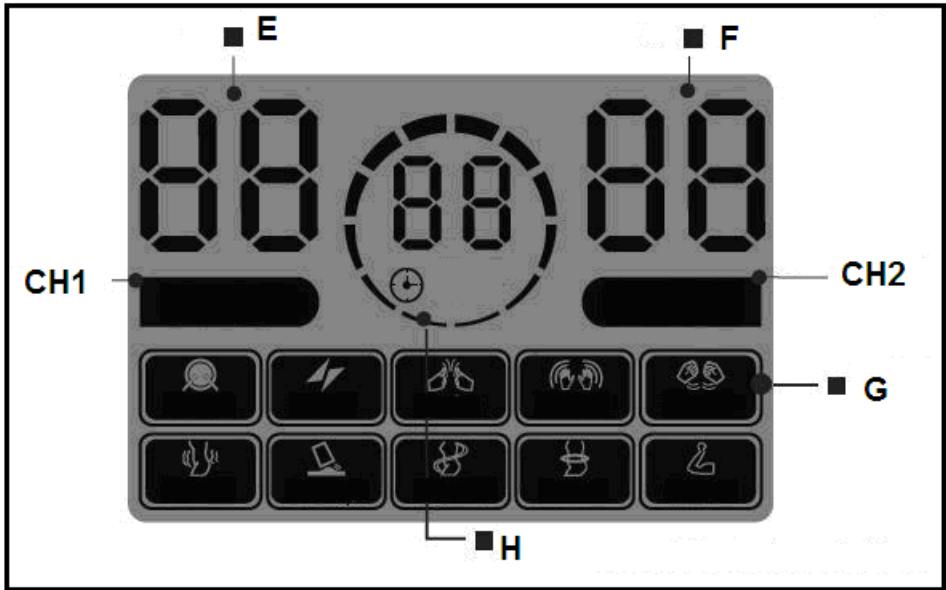E – Channel 1 intensity (40 intensity levels)

F – Channel 2 intensity (40 intensity levels)

G - Mode

H - Time (10 minutes).


CH1 – Channel 1

CH2 – Channel 2


BEFORE USE

To fit the batteries, remove the cover at the back of the remote control unit and carefully insert 2 «AA» type batteries into the battery compartment, making sure that the batteries' polarity is correct (matching the battery's «+» terminal with the «+» sign inside the battery compartment, and the «-» terminal with the «-» sign). Refit the battery cover.

1. Apply plenty of electro-stimulation gel on the electrodes.
2. Put the belt around your waist (the electrodes have to be in contact with skin. Otherwise there will be no stimulation). Adjust the width of the belt.

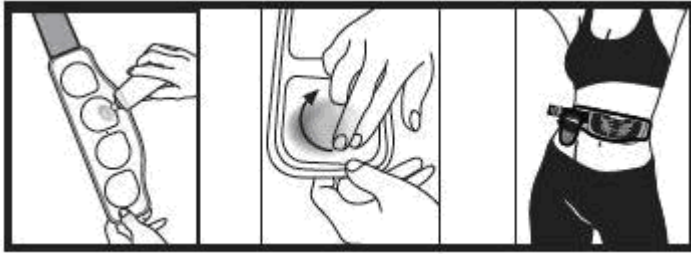3. Press the On/Off/ button (5) (Fig. 1)

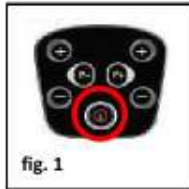4. The digital screen will switch on. Select the mode by pressing the mode selection button P+ or P- (6) and (7) (fig. 2)

5. Set the intensity for channel 1 by pressing the + or – buttons (8). Set the intensity for channel 2 by pressing the + or – buttons (9). (fig. 3)

6. The timer will start to countdown. The appliance will switch off automatically after 10 minutes of use.
7. You can change the mode at any time by pressing the mode selection button P+ or P- (6) and (7) (fig. 2). By default, the intensity will switch back to 0. Adjust the intensity by pressing the + or - buttons (8) or (9). (fig. 3)
8. You can switch the appliance off at any time. To do this, press and hold the On/Off button (5).

NOTE: If the electrodes are not in contact with the skin, the intensity will switch to 0 and the appliance automatically switches off.

POSSIBLE PROBLEMS

PROBLEM	CAUSE	SOLUTION
THE SCREEN DOES NOT SWITCH ON	<ul style="list-style-type: none"> - THERE ARE NO BATTERIES INSTALLED - THE BATTERIES ARE NOT INSTALLED CORRECTLY - BATTERIES ARE DEAD 	<ul style="list-style-type: none"> - FIT BATTERIES - CHECK THE POLARITY - REPLACE THE BATTERIES.
THE SCREEN SWITCHES ON BUT STIMULATION IS AT 0	<ul style="list-style-type: none"> - THE ELECTRODES ARE NOT IN CONTACT WITH THE SKIN - INTENSITY LEVEL AT 0 - SOFTWARE FROZEN - LOW BATTERY POWER 	<ul style="list-style-type: none"> - APPLY GEL AND PLACE ELECTRODES IN CONTACT WITH THE SKIN - INCREASE THE INTENSITY - RESET THE APPLIANCE BY REMOVING THE BATTERIES REFIT THE BATTERIES AFTER 1 MINUTE - REPLACE THE BATTERIES.
THE SCREEN SWITCHES ON BUT STIMULATION IS NOT NOTICEABLE	<ul style="list-style-type: none"> - LOW INTENSITY - NOT ENOUGH GEL ON THE ELECTRODES - GEL DIFFERENT TO THE ONE THAT COMES WITH THE PRODUCT USED - THE BELT HAS NOT BEEN FITTED CORRECTLY - THE ELECTRODES ARE DIRTY - LOW BATTERY POWER 	<ul style="list-style-type: none"> - INCREASE THE INTENSITY - APPLY MORE GEL - APPLY THE GEL THAT COMES WITH THE PRODUCT - FIT THE BELT CORRECTLY - CLEAN THE ELECTRODES - REPLACE THE BATTERIES.
THE BUTTONS DO NOT FUNCTION CORRECTLY	<ul style="list-style-type: none"> - SOFTWARE FROZEN - LOW BATTERY POWER 	<ul style="list-style-type: none"> - RESET THE APPLIANCE BY REMOVING THE BATTERIES REFIT THE BATTERIES AFTER 1 MINUTE - REPLACE THE BATTERIES.

MAINTENANCE

1. Clean the electrodes after using a soft cloth lightly moistened with water.
2. Make sure that no water gets inside the appliance.
3. Remove the batteries from the unit when it will not be used for prolonged periods. Leaking batteries may damage the unit.
4. Used batteries must be deposited in special collection containers for toxic waste. The law requires proper disposal of batteries.
5. Store the appliance in a dry, well ventilated place.

TECHNICAL SPECIFICATIONS

MODEL	BT BELT ELECTROSTIMULATOR YR30
FREQUENCY	1Hz – 60Hz
BATTERIES	2X1.5V AA ALKALINE BATTERIES
DIMENSIONS ELECTRODE SECTION	410X135MM
EXTENDABLE BELT	700X38MM
WEIGHT	300g

Warranty

This product is guaranteed for 2 years in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25th May, 1999, regarding certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10th July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT n° B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product.

Date of purchase.-	Purchased in.-
Purchaser.-	Dealer stamp.-

Please keep the receipt as proof of purchase.

Should you have any queries please do not hesitate to phone our customer service line:

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BH China Co., Ltd.